

For Research Use

TAKARA

**PrimeScript™ RT Master Mix
(Perfect Real Time)**

Product Manual

Table of Contents

I.	Description.....	3
II.	Components	3
III.	Materials Required but not Provided	4
IV.	Storage	4
V.	Features.....	4
VI.	Precautions for Use.....	4
VII.	Protocol.....	5
VIII.	Real-Time PCR.....	6
IX.	Experimental Example	9
X.	Appendix	10
XI.	Related Products	12

I. Description

This product is a reverse transcription reagent kit designed to perform reverse transcription optimized for 2 step real-time RT-PCR (RT-qPCR). It contains a 5X pre-mixed reagent containing all of the components needed for quantitative RT-PCR reverse transcription (PrimeScript RTase, RNase Inhibitor, Random 6 mers, Oligo dT Primer, dNTP Mixture, and reaction buffer), and a reaction can be started simply by adding template RNA and water. Because it uses PrimeScript RTase, which features excellent extendibility, this product, like PrimeScript RT reagent Kit (Perfect Real Time) (Cat. #RR037A/B), makes it possible to synthesize template cDNA for real-time PCR efficiently in a short time.

The cDNA obtained with this product can be used in both intercalator qPCR assay and probe qPCR assay. Please use this product in combination with quantitative PCR reagent such as TB Green® *Premix Ex Taq*™ II (Tli RNaseH Plus) (Cat. #RR820A/B)*¹, TB Green Fast qPCR Mix (Cat. #RR430A/B)*¹, and Probe qPCR Mix (Cat. #RR391A/B)*², depending upon your objectives.

*¹ We have begun the process of changing the names for Takara Bio's intercalator-based real-time PCR (qPCR) products to the "TB Green series". These products can be used the same way as before, as only the names are changing. Catalog number and product performance are unaffected by this transition.

*² Not available in all geographic locations. Check for availability in your area.

II. Components (200 reactions, 10 μl per reactions)

1. 5X PrimeScript RT Master Mix (Perfect Real Time)* ³	400 μl
2. RNase Free dH ₂ O	1 ml x 2
3. EASY Dilution (for Real Time PCR)* ⁴	1 ml

*³ Contains PrimeScript RTase, RNase Inhibitor, Oligo dT Primer, Random 6 mers, dNTP Mixture, and reaction buffer (containing Mg²⁺).

*⁴ To be used when preparing serial dilutions of total RNA and cDNA. EASY Dilution (for Real Time PCR) makes it possible to obtain a precise dilution down to very low concentrations. Moreover, this solution will not affect reverse-transcription and PCR reactivity. You can use the diluted solution directly in reverse-transcription and PCR reactions.

EASY Dilution (for Real Time PCR) may also be purchased separately (Cat. #9160).

Note: Please use EASY Dilution in combination with Takara Bio's real-time PCR Reagents. Compatibility with products from other manufacturers has not yet been verified.

III. Materials Required but not Provided

- Thermal cycler (or 37°C water bath and 85°C heat block)
- 0.2 ml and 1.5 ml microtubes (for reverse-transcription reaction)
- Micropipettes and tips (autoclaved)

IV. Storage -20°C

V. Features

1. The reaction can be started simply by adding template RNA and water.
2. Template cDNA for real-time PCR can be efficiently synthesized in a short time. This kit is best suited for 2 step real-time RT-PCR.
3. Optimal template cDNA for real-time PCR can be synthesized with 2 kinds of reverse-transcription primers: Oligo dT Primer and Random 6 mer.
4. A standard curve must be generated for the quantification of real-time RT-PCR. To generate a precise standard curve, it is important to dilute the total RNA and reverse-transcribed cDNA precisely to low concentrations. However, dilution of templates with water or TE buffer may narrow the range of the curve due to template instability at low concentrations. Using EASY Dilution (for Real Time PCR) makes for more accurate results at low concentrations, facilitating the creation of a wide-range standard curve.

VI. Precautions for Use

Read these precautions before use and follow them when using this product.

1. Briefly centrifuge the 5X PrimeScript RT Master Mix before use and allow the reagent to settle to the bottom of the tube. Since the 5X PrimeScript RT Master Mix has high viscosity, you will need to perform pipetting slowly and carefully, and also mix the reaction solution well.
2. When dividing up the reagent, please be sure to use new disposable tips to avoid contamination among samples.
3. This product is a completely pre-mixed reagent with reverse-transcription primers already added (Oligo dT Primer and Random 6 mers).

Note: This product cannot be used for reverse-transcription reactions using Gene-Specific Primer.

VII. Protocol**Reverse transcription reaction**

For RNA preparation method, please see Section X.

1. Prepare the reverse-transcription reaction solution on ice.

<Per reaction>

Reagent	Volume	Final conc.
5X PrimeScript RT Master Mix (Perfect Real Time)	2 μ l	1X
total RNA*		
RNase Free dH ₂ O	up to 10 μ l	

- * The scale of the reverse-transcription reaction can be increased as necessary. Reverse transcription of as much as 500 ng of total RNA is possible with 10 μ l of reaction solution.

2. Perform the reverse-transcription reaction after gently mixing the reaction solution.

37°C 15 min (reverse-transcription)

85°C 5 sec (for heat inactivation of reverse transcriptase)

4°C

Note: When the reverse-transcription solution obtained in Step 2 is to be introduced to a real-time PCR reaction, the volume of the solution should be 10% or less of the reaction volume for real-time PCR.

VIII. Real-Time PCR

The protocol for performing real-time PCR using TB Green *Premix Ex Taq II* (Tli RNaseH Plus) (Cat. #RR820A) after performing a reverse-transcription reaction with this kit is shown below.

[Protocol when using Applied Biosystems 7300/7500/7500 Fast Real-Time PCR System and StepOnePlus Real-Time PCR System]

* Please follow the procedures provided in the manual of the respective instrument.

1. Prepare the PCR mixture shown below.

<Per reaction>


Reagent	Volume	Volume	Final conc.
TB Green <i>Premix Ex Taq II</i> (2X)	10 μ l	25 μ l	1X
PCR Forward Primer (Tli RNaseH Plus) (10 μ M)	0.8 μ l	2 μ l	0.4 μ M*1
PCR Reverse Primer (10 μ M)	0.8 μ l	2 μ l	0.4 μ M*1
ROX Reference Dye (50X) or Dye II (50X)*2	0.4 μ l	1 μ l	1X
RT reaction solution (cDNA solution)*3	2 μ l	4 μ l	
Sterile purified water	6 μ l	16 μ l	
Total	20 μ l*4	50 μ l*4	

- *1 A final primer concentration of 0.4 μ M is most likely to yield a good result. Nevertheless, if there is an issue with reactivity, try to find an optimal concentration between 0.2 and 1.0 μ M.
- *2 The concentration for ROX Reference Dye II (50X) is lower than that for ROX Reference Dye (50X).
- Use ROX Reference Dye II (50X) when performing analyses with Applied Biosystems 7500 Real-Time PCR System or 7500 Fast Real-Time PCR System.
 - Use ROX Reference Dye (50X) when using StepOnePlus or 7300 Real-Time PCR System.
- *3 In a 20 μ l reaction volume, it is preferable to use a quantity of cDNA corresponding to 10 pg - 100 ng of total RNA template. In addition, the reverse-transcription reaction solution volume should correspond to less than 10% of the PCR reaction solution volume.
- *4 Prepare in accordance with the recommended volume for each instrument.

2. Initiate the reaction.

The recommended protocol for PCR is the shuttle PCR standard protocol described below. Try this protocol first and optimize PCR conditions as necessary. Perform a 3 step PCR when using a primer with low T_m value or when a shuttle PCR is not feasible.

1) Applied Biosystems 7300/7500 Real-Time PCR System, StepOnePlus

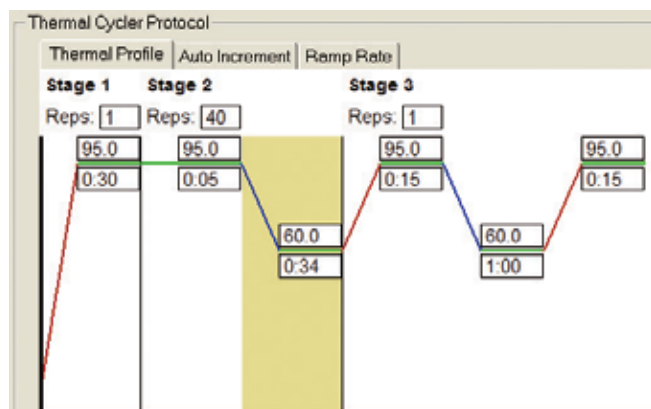

Shuttle PCR standard protocol

Stage 1: Initial denaturation

Reps: 1
95°C 30 sec

Stage 2: PCR
Reps: 40
95°C 5 sec
60°C 30 - 34 sec*

Dissociation stage

* With StepOnePlus,
set to 30 sec;
with 7300, set to 31 sec; and
with 7500, set to 34 sec

2) Applied Biosystems 7500 Fast Real-Time PCR System

Shuttle PCR standard protocol

Holding Stage

Number of cycle: 1
95°C 30 sec

Cycling Stage

Number of Cycles: 40
95°C 3 sec
60°C 30 sec

Melt Curve Stage

Note:

TaKaRa Ex Taq® HS used in this product is a hot-start PCR enzyme utilizing an anti-*Taq* antibody that suppresses polymerase activity. The initial denaturation step prior to PCR should be at 95°C for 30 sec. Longer heat treatment may decrease enzyme activity and affect amplification efficiency and quantification.

3. After the reaction is complete, check the amplification and melting curves and plot a standard curve if an assay will be performed.

Please refer to the instruction manual for your real time PCR instrument to read about analytical methods.

[Method Using Thermal Cycler Dice™ Real Time System //]

1. Prepare the PCR reaction mixture.

<Per reaction>


Reagent	Volume	Final conc.
TB Green <i>Premix Ex Taq</i> II (Tli RNaseH Plus) (2X)	12.5 μ l	1X
PCR Forward Primer (10 μ M)	1 μ l	0.4 μ M*1
PCR Reverse Primer (10 μ M)	1 μ l	0.4 μ M*1
RT reaction solution (cDNA solution)	2 μ l*2	
Sterile purified water	8.5 μ l	
Total	25 μl	

*1 Good results are mostly obtained with a final primer concentration of 0.4 μ M, but when there is a problem with reactivity, it is best to consider an optimal concentration in the range of 0.2 - 1.0 μ M.

*2 It is preferable to use a quantity of cDNA corresponding to 10 pg - 100 ng of total RNA template. In addition, the reverse-transcription solution introduced should correspond to 10% or less of the PCR reaction solution volume.

2. Start the reaction.

Shuttle PCR standard protocol (below) is recommended. Try this protocol first, and optimize the reaction condition if needed. When the shuttle protocol is difficult due to a primer with low Tm value, etc., try a 3 step PCR protocol.

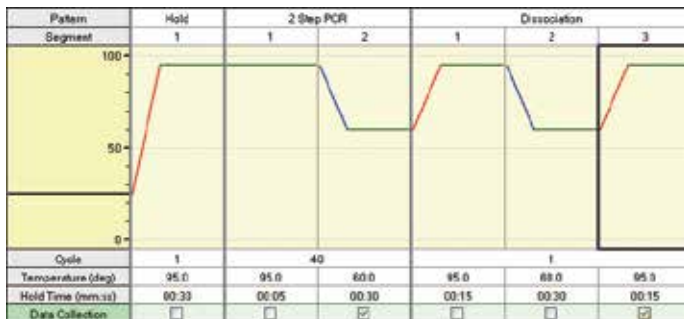

Shuttle PCR standard protocol

- Hold (initial denaturation)
 - Cycle: 1
 - 95°C 30 sec
- 2 Step PCR
 - Cycle: 40
 - 95°C 5 sec
 - 60°C 30 sec
- Dissociation

Note:
TAKARA Ex Taq HS used in this product is a hot-start PCR enzyme utilizing an anti-*Taq* antibody that suppresses polymerase activity. The initial denaturation step prior to PCR should be at 95°C for 30 sec. Longer heat treatment may decrease enzyme activity and affect amplification efficiency and quantification.

3. After the reaction is complete, check the amplification and melting curves and plot a standard curve if a quantitative assay will be performed.
When using Thermal Cycler Dice Real Time System //, please refer to its instruction manual to read analytical methods.

IX. Experimental Example


Reverse-transcription reaction time and amount of cDNA synthesized.

[Method]


1. Reverse-Transcription Reaction
 - Reagent: PrimeScript RT Master Mix (Perfect Real Time)
 - Template: Human Placenta total RNA 2 pg - 2 μg and sterilized water
 - Reaction volume: 20 μl
 - Reaction conditions: 37°C for 15, 30, or 60 min → 85°C for 5 sec → 4°C
2. Real Time PCR
 - Reagent: TB Green Premix Ex Taq II (Perfect Real Time)
 - Template: 2 μl of each of the aforementioned reverse-transcription reactions
 - Reaction volume: 25 μl
 - Target gene: *ACTB*
 - Reaction conditions: Standard protocol for use with Thermal Cycler Dice Real Time System

[Results]

Amplification curve

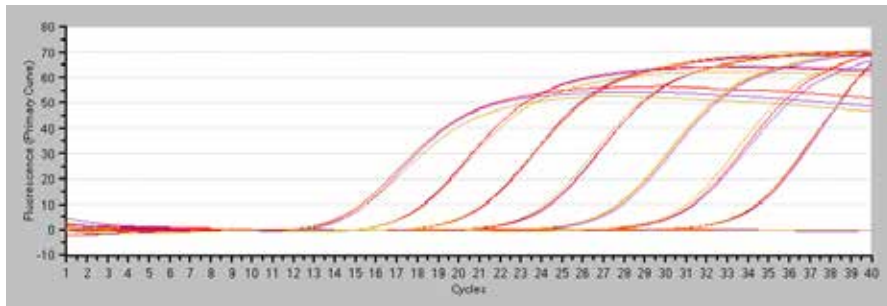

Standard curve

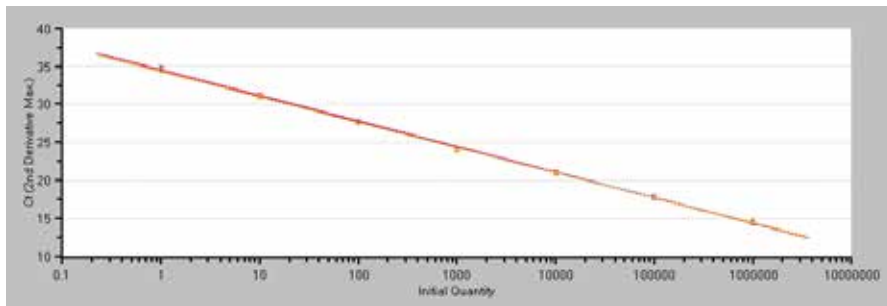

	Time	RSq	Eff (%)	Standard Curve
Violet	15 min	0.999	99.0	Y = -3.346 * LOG(X)+34.52
Red	30 min	1.000	98.9	Y = -3.349 * LOG(X)+34.46
Yellow	60 min	0.999	100.9	Y = -3.301 * LOG(X)+34.24

Reactivity was compared for cases where the reaction time was set at 15, 30, or 60 min. This experiment showed that reactions could proceed with equivalent efficiency over a wide range of template concentrations regardless of the reaction time.

X. Appendix

Preparation of RNA Samples

It is important to use a highly pure RNA samples for better cDNA yield. It is essential to inhibit RNase activity in the cells and also to prevent RNase derived from equipment and solutions used. Extra precautions should be taken during the sample preparation, including use of clean disposable gloves, dedication of a exclusively table for exclusive use for RNA preparation, and avoiding unnecessary speaking during assembly to prevent RNase contamination from operator sweat or saliva.

A. Equipment

Disposable plastic equipment should be used. For glass tools should be treated with the following procedure prior to use.

- (1) Hot-air sterilization (180°C, 60 min)
- (2) Sterilize glassware in 0.1% diethyl pyrocarbonate (DEPC) solution at 37°C for 12 hours. Autoclave (120°C, 30 min) to remove any residual DEPC.

It is recommended that all the equipment be used exclusively for RNA preparation.

B. Reagent

All reagents to be used in this experiment must be prepared using tools which were treated as described in previous section (Hot-air sterilization (180°C, 60 min) or DEPC treatment), and all purified water must be treated with 0.1% DEPC and autoclaved.

All reagents and purified water should be used exclusively for RNA experiments.

C. Method of Preparation for RNA Samples

Use of highly purified RNA obtained by GTC (Guanidine thiocyanate) method, etc. is recommended. RNA isolation kits such as RNAiso Plus (Cat. #9108/9109) can also be used for isolating high purity total RNA. The purified RNA sample should be dissolved in sterilized purified water or sterilized TE buffer.

D. Contamination with Genomic DNA and its countermeasure

In some cases, a total RNA sample may contain a small amount of genomic DNA, which could potentially become a PCR template. This could result in inaccurate results.

To avoid this to situation, either design primers which will not amplify genomic DNA or remove genomic DNA by DNase I treatment.

- (1) Designing a primer which will not amplify genomic DNA:

It is possible to design a primer that will not amplify genomic origin of DNA based on the exon and intron structure of genomic DNA. First, confirm the target genomic structure, and select a large intron region. Then design primers on the upstream and downstream sides of the intron region. If the intron is large enough, genomic DNA amplification will not occur. Even if size of the intron is not large, products resulting from amplification of genomic templates would be larger than products originating from cDNA, allowing identification by dissociation curve analysis. However, this approach cannot be applied for single-exon genes or genes that contain pseudogene. Moreover, species without introns or new species without well-identified genomic structures could demonstrate similar problems. In this case, we recommend performing the DNase I treatment described in (2).

- (2) After extraction of total RNA, remove genomic DNA by Recombinant DNase I (RNase-free) (Cat. #2270A) treatment. After the reaction, DNase I should be inactivated by either heat treatment or phenol/chloroform extraction.

Procedure

1. Prepare the following reagents:

Reagent	Volume
total RNA	20 - 50 μ g
10X DNase I Buffer	5 μ l
RNase Inhibitor	20 U
DNase I (RNase-free)	2 μ l (10 U)
DEPC treated water	up to 50 μ l

2. Incubate at 37°C for 20 min.
3. Perform one of the following procedures to inactivate DNase I.
 - A. Heat treatment
 - (1) Add 2.5 μ l of 0.5 M EDTA, and incubate at 80°C for 2 min.
 - (2) Increase reaction volume to 100 μ l with DEPC treated water.
 - B. Phenol/Chloroform extraction
 - (1) Mix 50 μ l of DEPC treated water and 100 μ l of phenol/chloroform/isoamyl alcohol (25 : 24 : 1) together.
 - (2) Centrifuge at 15,000 rpm for 5 min at room temperature then transfer the upper layer to a new tube.
 - (3) Add equal amount of Chloroform/isoamyl alcohol (24 : 1) and mix.
 - (4) Centrifuge at 15,000 rpm for 5 min at room temperature then transfer the upper layer to new tube.
4. Add 10 μ l of 3 M sodium acetate and 250 μ l of cold ethanol, then incubate on ice for 10 min.
5. Centrifuge at 15,000 rpm for 15 min at 4°C and remove the top clear liquid.
6. Wash the precipitate with 70% ethanol, and then centrifuge at 15,000 rpm for 5 min at 4°C and remove the supernatant.
7. Dry the precipitate.
8. Dissolve the precipitate in the proper amount of DEPC treated water.

E. Confirmation of Genomic DNA Contamination

The presence of genomic DNA can be verified by real-time PCR without reverse transcription. It is convenient to use a primer that can do PCR amplification from both genomic DNA and mRNA for this experiment, and one tube should include an exogenous genomic DNA aliquot. In addition, a primer which was designed not to amplify genomic DNA could amplify derived from pseudogene. In this case, this procedure can also be use to verify DNA.

XI. Related Products

TB Green® *Premix Ex Taq*™ II (Tli RNaseH Plus) (Cat. #RR820A/B)

TB Green® Fast qPCR Mix (Cat. #RR430A/B)*

TB Green® *Premix Ex Taq*™ (Tli RNaseH Plus) (Cat. #RR420A/B)

Probe qPCR Mix (Cat. #RR391A/B)*

EASY Dilution (for Real Time PCR) (Cat. #9160)

Thermal Cycler Dice™ Real Time System III (Cat. #TP950/TP970/TP980/TP990)*

Thermal Cycler Dice™ Real Time System II (Cat. #TP900/TP960)*

Thermal Cycler Dice™ Real Time System *Lite* (Cat. #TP700/TP760)*

PrimeScript™ RT Reagent Kit (Perfect Real Time) (Cat. #RR037A/B)

PrimeScript™ RT reagent Kit with gDNA Eraser (Perfect Real Time) (Cat. #RR047A/B)

RNAiso Plus (Cat. #9108/9109)

* Not available in all geographic locations. Check for availability in your area.

Takara Ex Taq and TB Green are registered trademarks of Takara Bio Inc.

PrimeScript, Thermal Cycler Dice, and *Premix Ex Taq* are trademarks of Takara Bio Inc.

NOTE: This product is for research use only. It is not intended for use in therapeutic or diagnostic procedures for humans or animals. Also, do not use this product as food, cosmetic, or household item, etc.

Takara products may not be resold or transferred, modified for resale or transfer, or used to manufacture commercial products without written approval from Takara Bio Inc.

If you require licenses for other use, please contact us by phone at +81 77 565 6973 or from our website at www.takara-bio.com.

Your use of this product is also subject to compliance with any applicable licensing requirements described on the product web page. It is your responsibility to review, understand and adhere to any restrictions imposed by such statements.

All trademarks are the property of their respective owners. Certain trademarks may not be registered in all jurisdictions.
